


# Easy HiFlype

## NITINOL SELF EXPANDING STENT FOR ILIAC VESSELS

Easy HiFlype is a dedicated self expanding stent with hybrid cell design maximizing radial force and precision for iliac vessels.


- Proven bio and haemo-compatibility enabled by the Bio Inducer Surface coating, accelerating strut coverage, reducing thrombogenicity and sealing against heavy metal ions release
- Uniform radial force and excellent scaffolding with a unique hybrid cell architecture
- Superior track and high precision upon stent release provided by a unique single-hand delivery system

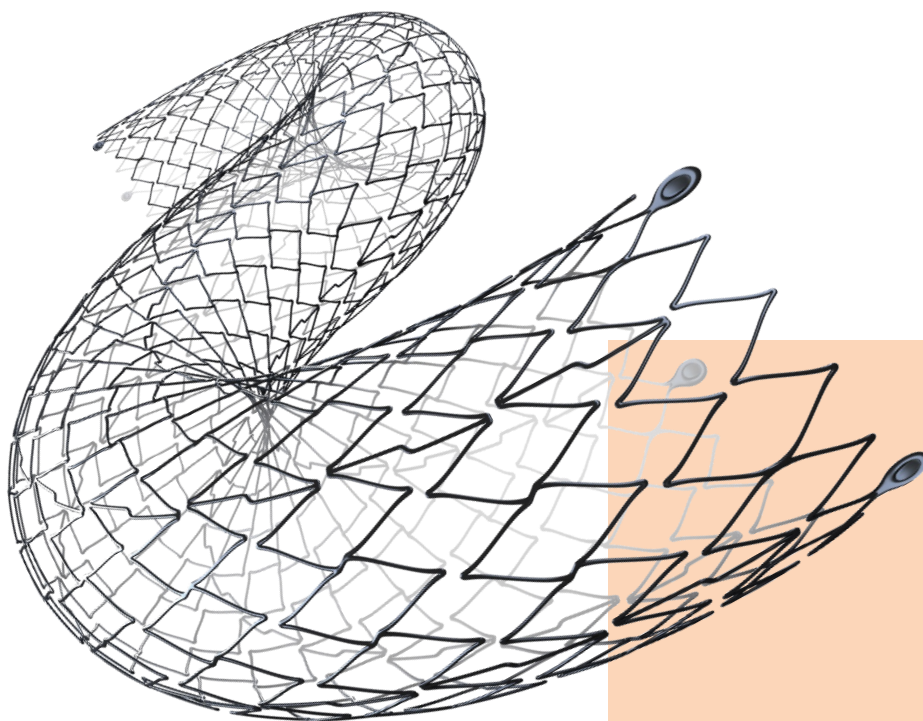


# Easy HiFlype

NITINOL SELF EXPANDING STENT FOR ILIAC VESSELS


## UNIQUE FEATURES

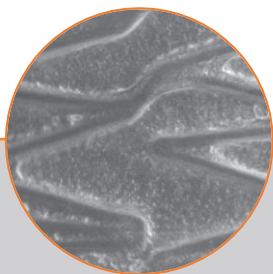

BIS coating accelerates the rate of endothelialization, reduces thrombogenicity and sealing against heavy metal ions release

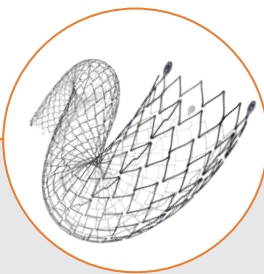


Uniform radial force and excellent scaffolding with a unique hybrid cell architecture

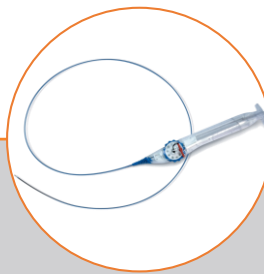

Superior track and high precision upon stent release provided by a unique single-hand delivery system

## TECHNICAL DATA

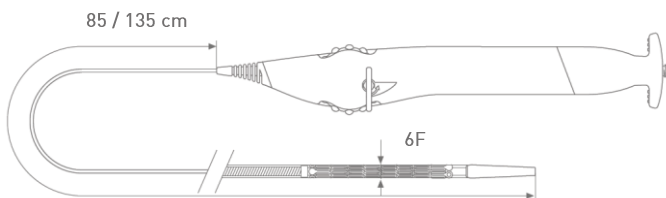

CATHETER SPECIFICATIONS		STENT SPECIFICATIONS	
Category:	Over The Wire (OTW)	Category:	Bio Inducer Surfaced Self Expanding BMS
Guiding catheter compatibility:	6F	Material:	Nickel Titanium Alloy (nitinol)
Guidewire compatibility:	0.035"	Strut Coating	Bio Inducer Surface
Usable catheter length:	85 - 135 cm	Stent diameters	9 - 12 mm
External Sheath:	Braided Polyamide	Stent lengths	20 - 100 mm
Dual Release System:	Thumbwheel for gradual deployment Pull back for rapid deployment	X-ray visibility:	Six dot tantalum markers

## ORDER INFORMATION EASY HIFLYPE

Usable Catheter Length	L [mm] ∅ [mm]	20	40	60	80	100
		ICEF9020S	ICEF9040S	ICEF9060S	ICEF9080S	ICEF90100S
85 cm	9	ICEF10020S	ICEF10040S	ICEF10060S	ICEF10080S	ICEF100100S
	10		ICEF12040S	ICEF12060S	ICEF12080S	ICEF120100S
	12					
135 cm	9	ICEF9020L	ICEF9040L	ICEF9060L	ICEF9080L	ICEF90100L
	10	ICEF10020L	ICEF10040L	ICEF10060L	ICEF10080L	ICEF100100L
	12		ICEF12040L	ICEF12060L	ICEF12080L	ICEF120100L

[www.alvimedica.com](http://www.alvimedica.com)

Please contact your local sales representative for product availability. The images shown within this document are of a purely illustrative nature and are not intended to depict reality.

Please consult product labels and package inserts for indications, contraindications, hazards, warnings, cautions and instruction for use.

©2020, Alvimedica. All Rights Reserved.

Manufactured by: CID S.p.A.  
member of Alvimedica group  
Via Crescentino, sn - 13040 Saluggia (VC), Italy

