


Easy Flype

NITINOL SELF EXPANDING STENT FOR SFA

Easy Flype is a dedicated self expanding stent with hybrid cell design maximizing flexibility and conformability in SFA.


- Proven bio and haemo-compatibility enabled by the Bio Inducer Surface coating, accelerating strut coverage, reducing thrombogenicity and sealing against heavy metal ions release
- Effective balance between flexibility (for optimal conformability) and radial force (for highly calcified vessels), available up to 150 mm in length
- Superior track and high precision upon stent release provided by a unique single-hand delivery system

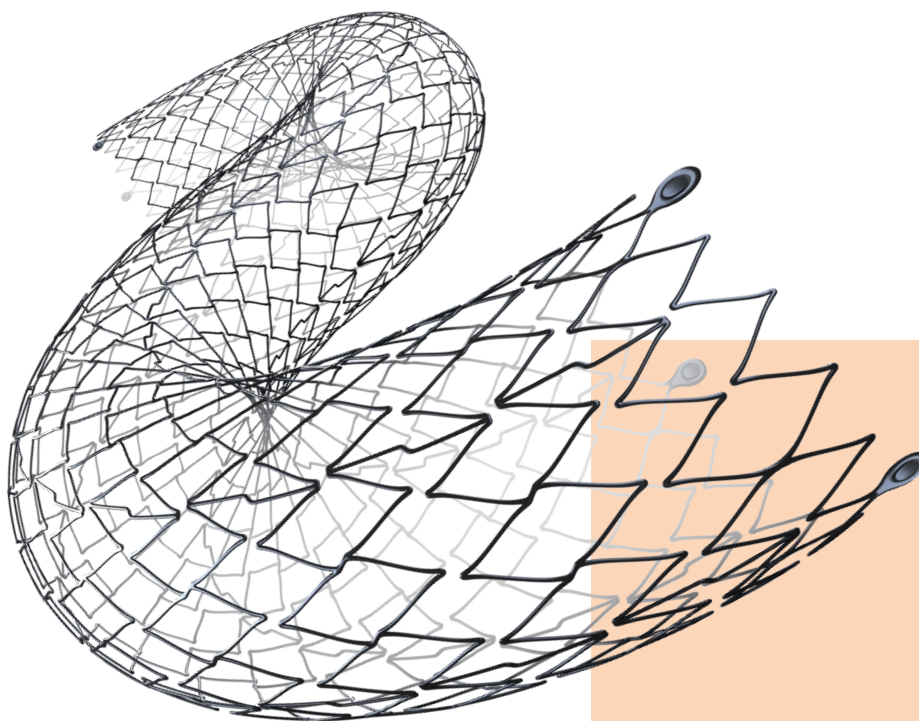


Easy Flype

NITINOL SELF EXPANDING STENT FOR SFA


UNIQUE FEATURES

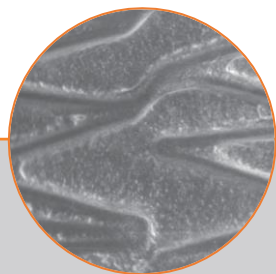

BIS coating accelerates the rate of endothelialization, reduces thrombogenicity and sealing against heavy metal ions release

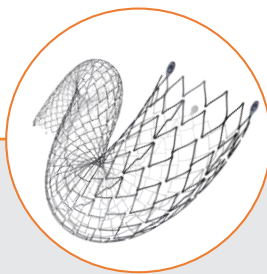


Effective balance between flexibility (for optimal conformability) and radial force (for highly calcified vessels), available up to 150 mm in length

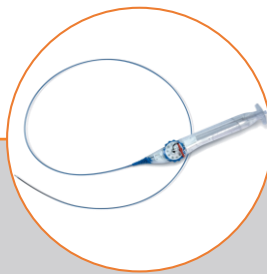

Superior track and high precision upon stent release provided by a unique single-hand delivery system

TECHNICAL DATA

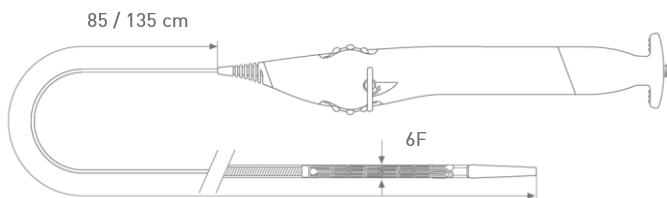

CATHETER SPECIFICATIONS		STENT SPECIFICATIONS	
Category:	Over The Wire (OTW)	Category:	Bio Inducer Surfaced Self Expanding BMS
Guiding catheter compatibility:	6F	Material:	Nickel Titanium Alloy (nitinol)
Guidewire compatibility:	0.035"	Strut Coating	Bio Inducer Surface
Usable catheter length:	85 - 135 cm	Stent diameters	6 – 8 mm
External Sheath:	Braided Polyamide	Stent lengths	20 – 150 mm
Dual Release System:	Thumbwheel for gradual deployment Pull back for rapid deployment	X-ray visibility:	Six dot tantalum markers

ORDER INFORMATION EASY FLYPE

Usable Catheter Length	L [mm] ∅ [mm]	20	40	60	80	100	120	150
		85 cm	ICEF6020S	ICEF6040S	ICEF6060S	ICEF6080S	ICEF60100S	ICEF60120S
	7	ICEF7020S	ICEF7040S	ICEF7060S	ICEF7080S	ICEF70100S	ICEF70120S	ICEF70150S
	8	ICEF8020S	ICEF8040S	ICEF8060S	ICEF8080S	ICEF80100S	ICEF80120S	ICEF80150S
135 cm	6	ICEF6020L	ICEF6040L	ICEF6060L	ICEF6080L	ICEF60100L	ICEF60120L	ICEF60150L
	7	ICEF7020L	ICEF7040L	ICEF7060L	ICEF7080L	ICEF70100L	ICEF70120L	ICEF70150L
	8	ICEF8020L	ICEF8040L	ICEF8060L	ICEF8080L	ICEF80100L	ICEF80120L	ICEF80150L

www.alvimedica.com

Please contact your local sales representative for product availability. The images shown within this document are of a purely illustrative nature and are not intended to depict reality.

Please consult product labels and package inserts for indications, contraindications, hazards, warnings, cautions and instruction for use.

©2020, Alvimedica. All Rights Reserved.

Manufactured by: CID S.p.A.
member of Alvimedica group
Via Crescentino, sn - 13040 Saluggia (VC), Italy

